

For collection, transport, maintenance
and long term freeze storage



**Advanced
Flocked Swabs**



Virus, Chlamydia, Mycoplasma and Ureaplasma

UTM

UNIVERSAL TRANSPORT MEDIUM

• **Room Temperature Stable** •

One transport - multiple tests - no sample rejection

- Culture - all viruses, Chlamydia, mycoplasma and ureaplasma
- Enzyme Immuno Assays (EIA)
- PCR
- Nucleic acid amplification assays including *Chlamydia trachomatis* and *Neisseria gonorrhoeae* (Ct/Ng).

Product stable at room temperature for 12 months

Organism viability maintained at room temperature and on ice

Choice of tube size, medium fill volume and product format


Variety of patient specimen collection kits available comprising tubes of media and swabs or bulk packs of media tubes only.

Advanced nylon flocked swabs increase sample release

UTM is supplied with unique flocked swabs which increase sample release and sensitivity in tissue culture or assays.


Precision snap molded breakpoint on swab shafts

Plastic swab applicators supplied with UTM kits feature exclusive molded breakpoints.

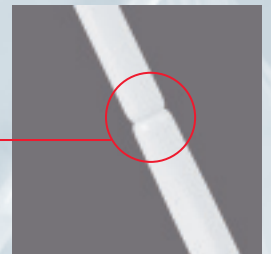

Capture cap docks and locks-in the swab stick


UTM tube caps incorporate a unique internal design which functions like a docking system automatically capturing the end of the swab stick and securing it firmly into the cap.

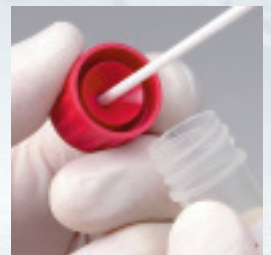

Internal conical shaped tube

Tubes stand upright on the bench without the need for tube racks.

Glass beads in each tube assist the release of sample from the swab and mixing prior to analysis when used with a vortex mixer.

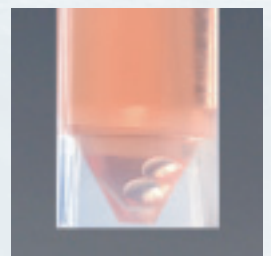

Copan Italia S.p.A.

Via F. Perotti, 10 • 25125 Brescia, Italy

Tel: (+39) 030 268 7211 • Fax: (+39) 030 268 7250

E-mail: customer-care@copanitalia.com web: www.copaninnovation.com

Copan UTM Ordering Information



Catalog No.	COPAN UTM - BULK TUBES (NO SWABS)	Packaging
330C	16x100 mm tube, 3 ml UTM medium (bulk pack)	6 x 50 kits
350C	12x80 mm mini tube, 1 ml UTM medium (bulk pack)	6 x 50 kits

Catalog No.	COPAN UTM - TUBE & FLOCKED SWAB	Packaging
346C	16x100 mm tube, 3 ml UTM medium, 1 Regular FLOCKED Swab, Molded Breakpoint	6 x 50 kits
307C	16x100 mm tube, 3 ml UTM medium, 1 Minitip FLOCKED Swab, Molded Breakpoint	6 x 50 kits
305C	16x100 mm tube, 3 ml UTM medium, 1 Flexible Nasopharyngeal FLOCKED Swab, Molded Breakpoint	6 x 50 kits
359C	12x80 mm mini tube, 1 ml UTM medium, 1 Regular FLOCKED Swab, Molded Breakpoint	6 x 50 kits
361C	12x80 mm mini tube, 1 ml UTM medium, 1 Minitip FLOCKED Swab, Molded Breakpoint	6 x 50 kits
360C	12x80 mm mini tube, 1 ml UTM medium, 1 Flexible Nasopharyngeal FLOCKED Swab, Molded Breakpoint	6 x 50 kits

Catalog No.	COPAN UTM - TUBE & STANDARD SWABS	Packaging
302C	16x100 mm tube, 3 ml UTM medium, 2 regular size plastic/polyester applicators	6 x 50 kits
303C	16x100 mm tube, 3 ml UTM medium, 2 plastic/polyester applicators (1 regular, 1 minitip)	6 x 50 kits
328C	16x100 mm tube, 3 ml UTM medium, 1 regular size plastic/polyester applicator	6 x 50 kits
335C	16x100 mm tube, 3 ml UTM medium, 3 plastic/polyester applicators (2 regular, 1 minitip)	6 x 50 kits
343C	16x100 mm tube, 1,5 ml UTM medium, no glass beads, 2 regular size plastic/polyester applicators	6 x 50 kits

Collection of Virus, Chlamydia, Mycoplasma and Ureaplasma

*For collection, transport,
maintenance and long
term freeze storage*



K37 06.07



Copan Italia S.p.a. Via F. Perotti, 10, Brescia, 25125, Italy

Tel: +39 030 268 7211 ¥ Fax: +39 030 268 7250

E-mail: info@copanitalia.com Check out our website: www.copaninnovation.com

Copan Flocked Swabs:
Worldwide Patent No. WO 2004/086979

